
N.168

# Plattenkondensator

## Parallel plate capacitor

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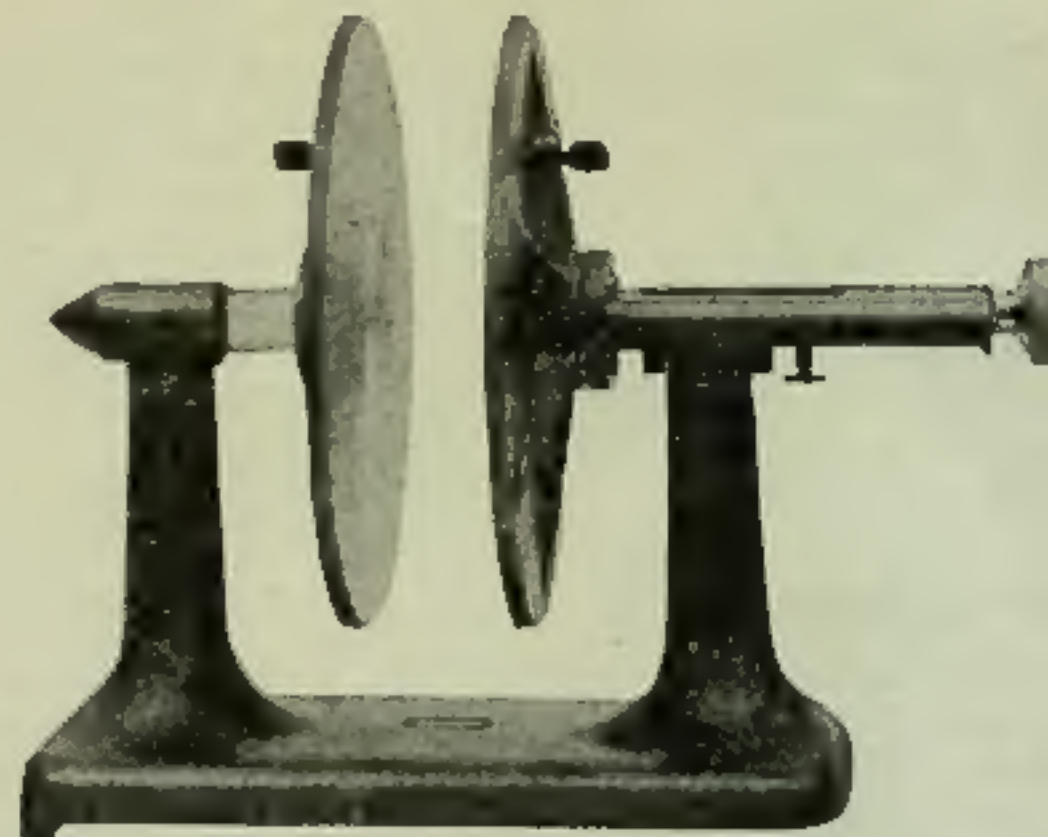The parallel plate capacitor is used in all experiments in the homogeneous electric field, that is, for introducing such quantities as quantity of charge, field strength, capacitance, charge density, and the phenomena associated with electrostatic induction. The permittivities of some substances and the dynamic effect in the electric field can be determined.

#### 1. Description

The parallel plate capacitor consists of a light metal frame supporting a stationary and a movable plate of 26 cm dia. and approx. 7 mm thickness each. It has two alternative three-point supports so that the plates can be used horizontally or vertically.

The plates are also made of light metal. A high-grade insulator separates the stationary plate from the frame. The movable plate can be

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#### Notes

- The five-digit numbers quoted in parentheses refer to the catalogue numbers of the respective apparatus.
- The data DC ... refer to the experiment descriptions published as a collection of leaflets in the "LEYBOLD PHYSICS LEAFLETS".  
Published to date:

	Cat. No.
1st series (72 sheets) . . . . .	599 21
2nd series (52 sheets) . . . . .	599 22
3rd series (49 sheets) . . . . .	599 23
4th series (50 sheets) . . . . .	599 24

- The specifications and illustrations are not binding in every detail for the design of the apparatus. It is our policy always to keep our manufacturing programme right up to date so that it makes full allowance for the latest knowledge acquired in all scientific and technical fields.

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adjusted by means of a covered post guide that can be arrested by a setscrew. The maximum plate spacing is 70 mm. For coarse adjustment of the plate spacing the setscrew should be loosened and the movable plate displaced by hand. The fine adjustment device controls a range of approx. 20 mm and can be operated by means of a small hand wheel when the setscrew is tight. The plate spacing can be read on a millimetre scale with vernier.

## 2. Operation

Additional equipment required for the measurements are an electrometer, e.g. electrostatic voltmeter for 1 kV (540 38), and the measuring amplifier (532 01) in combination with a demonstration multirange meter (531 86 or 87).

### a) Adjustment

Should it become necessary to correct the parallel adjustment of the plates, this can be done by means of the ball joint connecting the movable plate with the guide. In that case the three cheese-head screws around the ball joint should be loosened, and the movable plate should be pushed against the stationary one. Then the screws should be tightened a little, and the plates should be moved some distance apart and brought to contact again by means of the fine adjustment device. Should the plates then fail to make contact over the whole of their area, the screw next to the point of contact must be drawn tighter. Repeat this procedure until the plates are sufficiently parallel. Check by viewing the air gap between the two plates from several sides.

### b) Corrections

The values for the capacitance of the parallel plate capacitor found from the equation  $C_{th} = \epsilon_0 \frac{F}{d}$  ( $F$  = surface area of capacitor plate;  $d$  = plate spacing) must be increased by a capacitance  $C_R$  (from 7 up to 15 pF) due to the fringing of the field at the boundaries of the capacitor plates, and by a constant value  $C_K$  which is mainly due to the fringe effect towards the support, and amounts to approx. 7 pF.

Table 1 shows the capacitances and these corrections as a function of the plate spacing  $d$ .

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**c) Simple worked example**

The p. d. across a charged, insulated capacitor varies with the plate spacing. If the fringe effect at the boundaries and towards near objects is considered, it can be shown that  $U \sim d$  or  $\frac{U}{d} = \text{constant}$ .

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The p.d. should be measured by means of electrostatic instruments, the electrostatic voltmeter for 1 kV (540 38) being particularly suitable for the purpose on account of its very good insulation. In addition to the capacitances already listed, this arrangement involves further additional capacitances: the capacitance of the voltmeter ( $C_V$ ), which is 24 pF on an average, and the capacitance of the lead towards the voltmeter ( $C_S$ ), a connecting lead of 50 cm length generally having a total capacitance of approx. 4 pF.

Thus the total capacitance  $C$  of the circuit is composed of

$$C_{th} + C_R + C_K + C_V + C_S$$

where  $C_K$  and  $C_S$  are constant,  $C_V$  is only slightly variable, and even  $C_R$  can be considered nearly constant over small regions (voltage ratio 1:4).

The measurements plotted in diagram 3 show that there is a linear relation between  $\frac{1}{U}$  and  $\frac{1}{d}$ .

The slope  $a$  of the straight line depends on the geometry of the arrangement.

The correct values are obtained by displacing the coordinate system by a distance  $b$ :

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$$\frac{1}{U} - b = \frac{1}{U'}$$

Tableau 2

When using the corrected p.d.'s  $U'$ , the relation asserted at the beginning can be verified.

In order to improve the accuracy of measurement, the arrangement was initially charged by a source of D.C. which a previous check with a current-carrying voltmeter had shown to yield 300 V, the accuracy of reading of the high-tension voltmeter being insufficient in this range.

### 3. Notice

The high-grade insulation of the stationary plate makes it possible to experiment at leisure without risking substantial leakage losses. However, it is essential to keep the insulator

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clean, and Green's lens tissue (305 00) is particularly suitable for cleaning it.

When experimenting in humid air it is recommended to blow a hot air current against the capacitor.

At very small plate spacings (some tenths of a millimetre) and high p. d.'s the insulation may deteriorate on account of the plates being bridged by particles of dust and fibres that have settled between them. This fault can be overcome if both plates are covered with a very thin film of grease, or if the dust particles are burnt off with a Bunsen burner. But such small plate spacings will anyway only be used if a mirror galvanometer is the only available charge-measuring instrument.

#### 4. Experiments

The capacitance of a parallel plate capacitor with variable plate spacing.  
DC 537.224; a: 537.721.8, 4th series, leaflet 18

Capacitance determinations on the parallel plate capacitor  
Booklet: "The measuring amplifier and its application"

Experiments with a parallel plate capacitor  
Booklet: Wulf, "Das Wulf-Elektroskop"

Influence of the insulator material  
Booklet: Wulf, "Das Wulf-Elektroskop"

Force equals charge times electric field  
DC 537.211; a, 4th series, leaflet 14

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